


## INTRODUCTION


- The Comet Assay (Single Cell Gel Electrophoresis) is one of the most used techniques for measuring DNA damage<sup>1</sup> due to its simplicity, sensitivity and versatility<sup>2</sup>.
- The theoretical scientific basis of this methodology consists on the increased electrophoretic mobility of damaged DNA.
- It has a broad range of applications such as: compounds' genotoxicity/protective testing, DNA repair studies and ecological monitoring.
- With a huge amount of optimized protocols for superior eukaryotes, the Comet Assay has been poorly utilized among microorganisms.
- However, Azevedo et al. (2011)<sup>3</sup> developed an optimized protocol for *Saccharomyces cerevisiae* that reinforce the utility and potential of this method for microbial systems.
- Therefore, in this work we created a reproducible and optimized Comet Assay protocol for use in *Ashbya gossypii* by adaptation of the Yeast Comet Assay<sup>3</sup> (Fig. 1). This protocol allows the measuring of oxidative and non-oxidative DNA damage.**
- Moreover, as riboflavin overproduction in *A. gossypii* has been associated with stress (mainly oxidative), we went to check if this trait was also associated with increased DNA damage. Therefore, we used this recently created protocol to investigate the DNA damage accumulation in riboflavin overproducing and non-overproducing *A. gossypii* wild strains.**

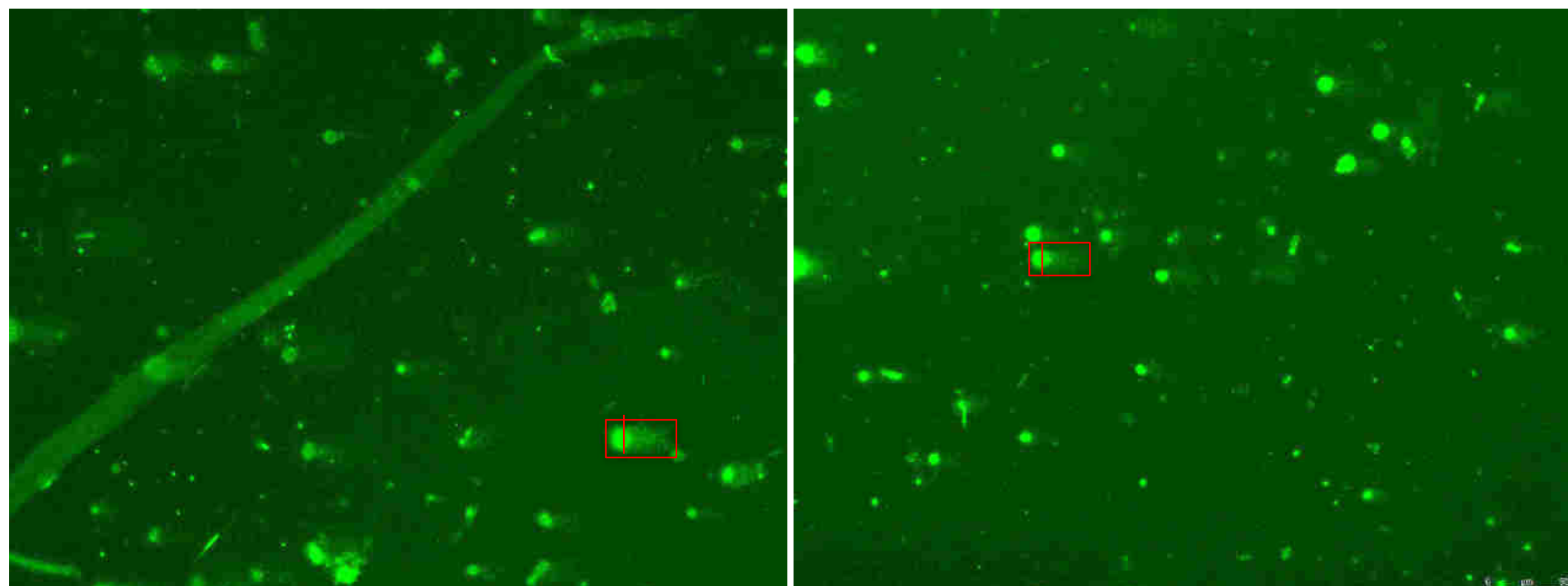

**Figure 1** - Representative fluorescent microscopy images resultant from the application of the *A. gossypii* Comet Assay developed in this work in cells treated with H<sub>2</sub>O<sub>2</sub>. The images were acquired at 200x magnification. White bar on the right image = 50 µm.

## METHODS


### 1. RADIAL GROWTH OF *A. gossypii* WILD STRAINS

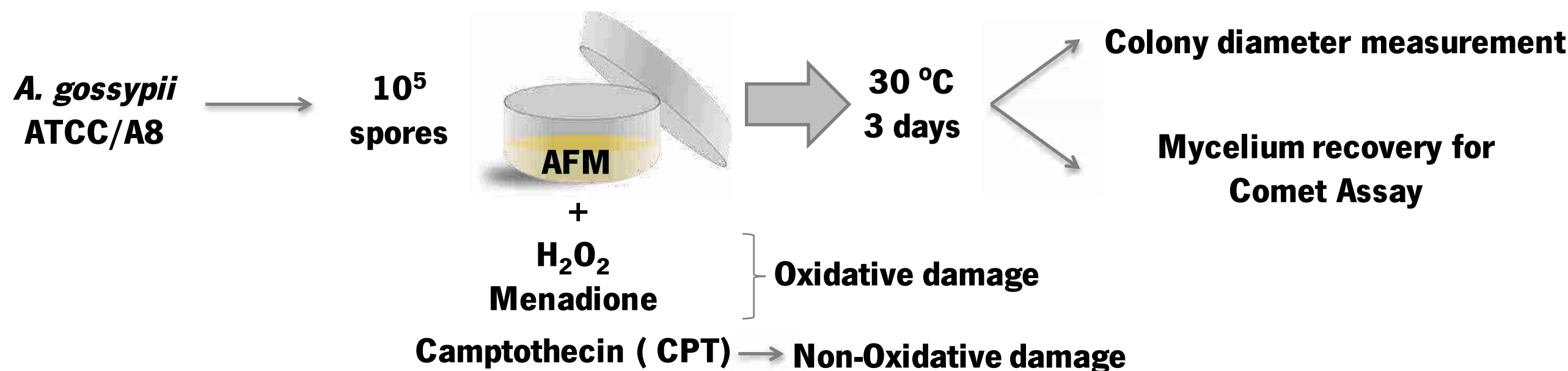


- A. gossypii* ATCC10895 (kindly provided by Prof. P. Philippsen, University of Basel) and *A. gossypii* A8 (kindly provided by Prof. J.L. Revuelta, University of Salamanca).
- For growth under light exposure a fluorescent lamp (LUMILUX T5 HE, OSRAM) was used with an irradiance level of 16 µmol m<sup>-2</sup> s<sup>-1</sup>, measured by a LI-250 Light Meter with a LI-190 quantum sensor (LI-COR, USA).

### 2. THE *A. gossypii* COMET ASSAY


#### 2.1 PROTOPLAST GENERATION

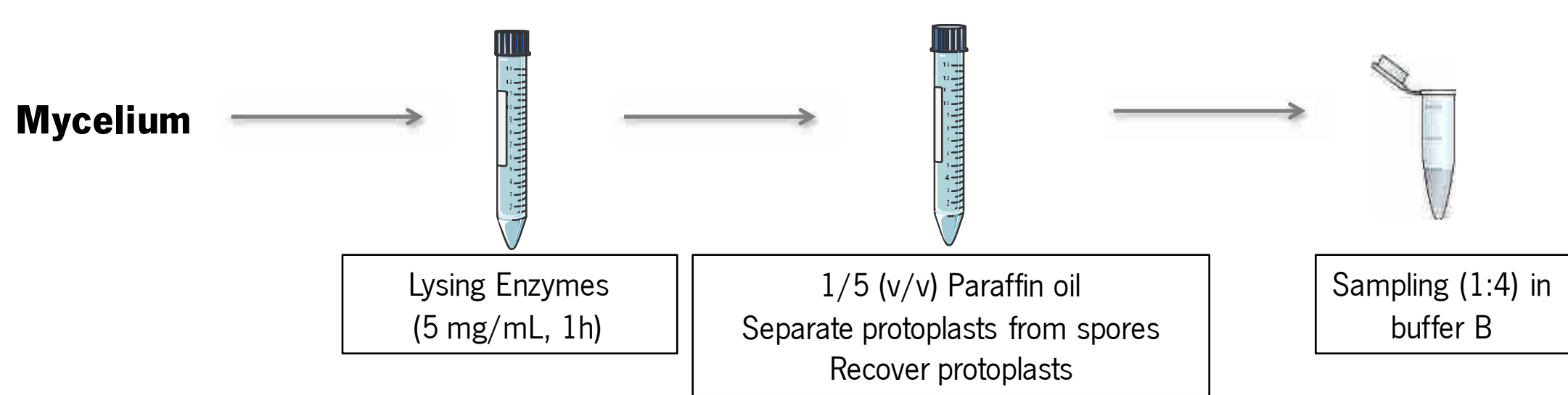

#### 2.2 BUFFER INCUBATIONS AND ELECTROPHORESIS

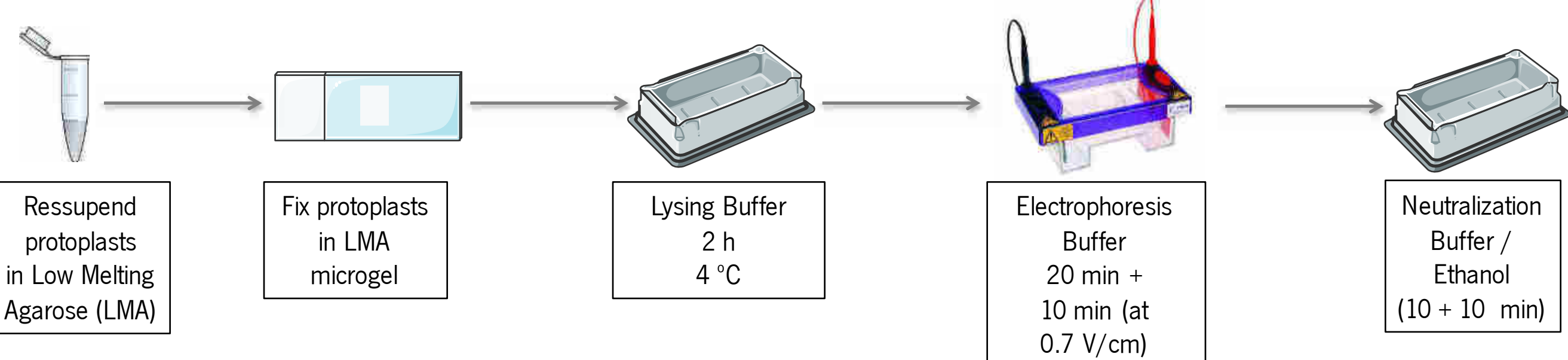


#### 2.3 COMETS VISUALIZATION AND ANALYSIS

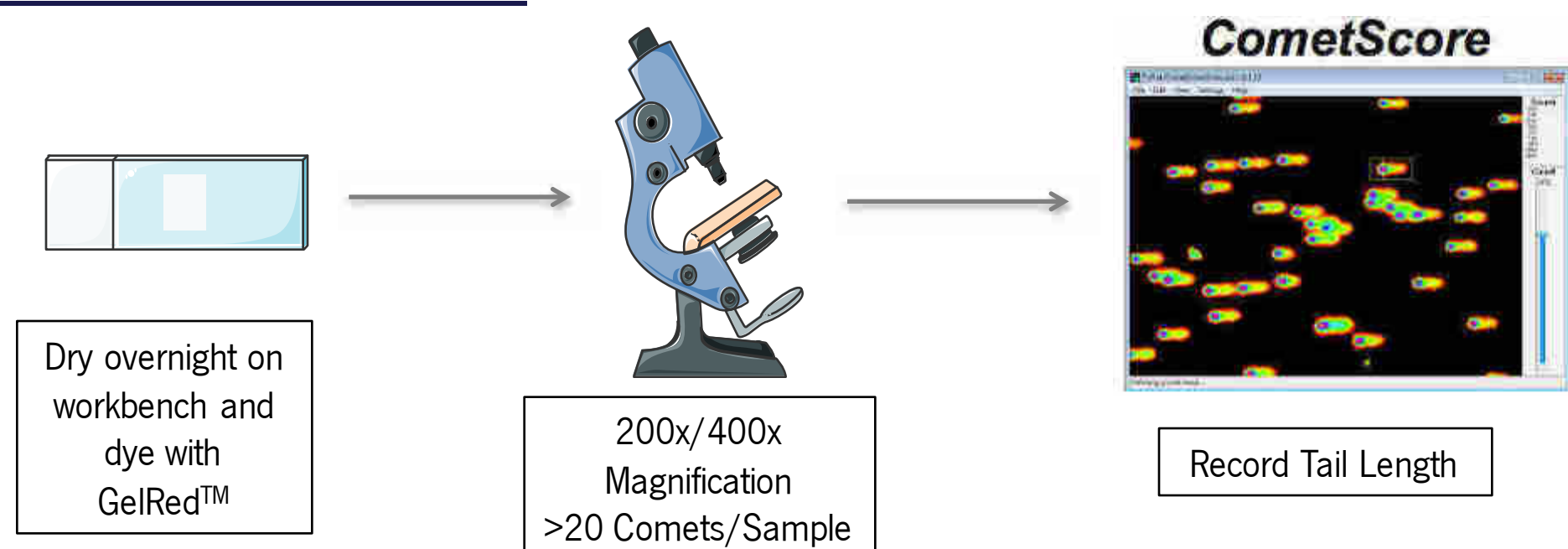


## RESULTS


### 1. THE *A. gossypii* COMET ASSAY IS ABLE TO MEASURE OXIDATIVE (H<sub>2</sub>O<sub>2</sub> AND MENADIONE) AND NON-OXIDATIVE (CPT) DNA DAMAGE

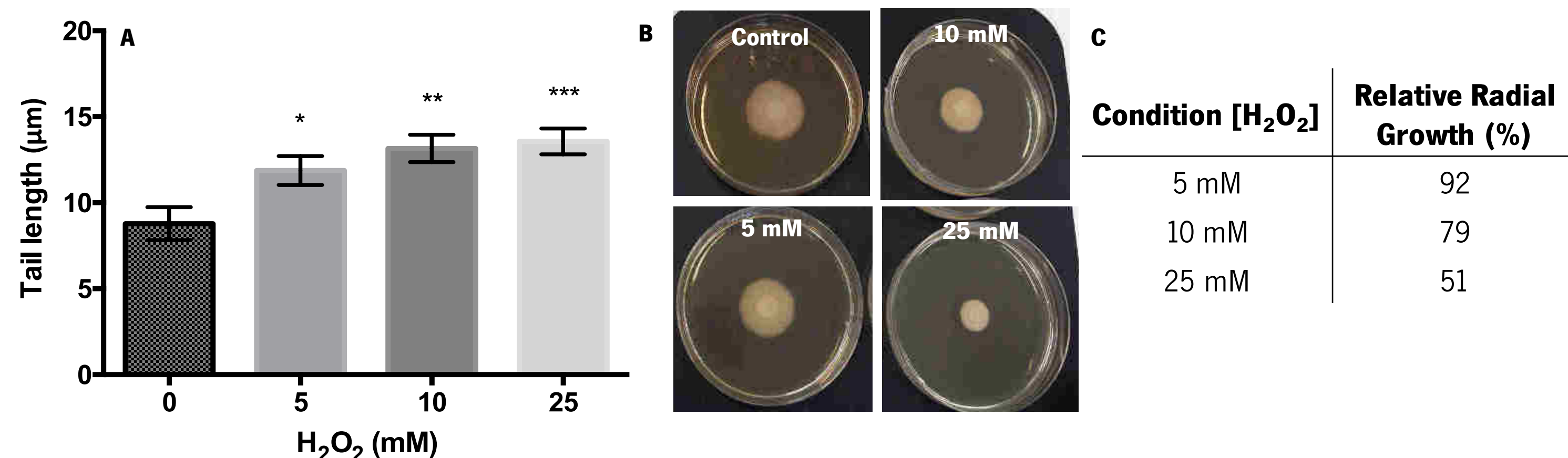

**Figure 2.1** – DNA damage and radial growth induced by exposure to H<sub>2</sub>O<sub>2</sub>. **(A)** DNA damage is represented as mean ± SEM of the tail length. **(B)** Representative photos of mycelia after growth for 3 days exposed to different concentrations of H<sub>2</sub>O<sub>2</sub>. **(C)** Relative growth of cells to H<sub>2</sub>O<sub>2</sub> in comparison to the control condition (0 mM H<sub>2</sub>O<sub>2</sub>).

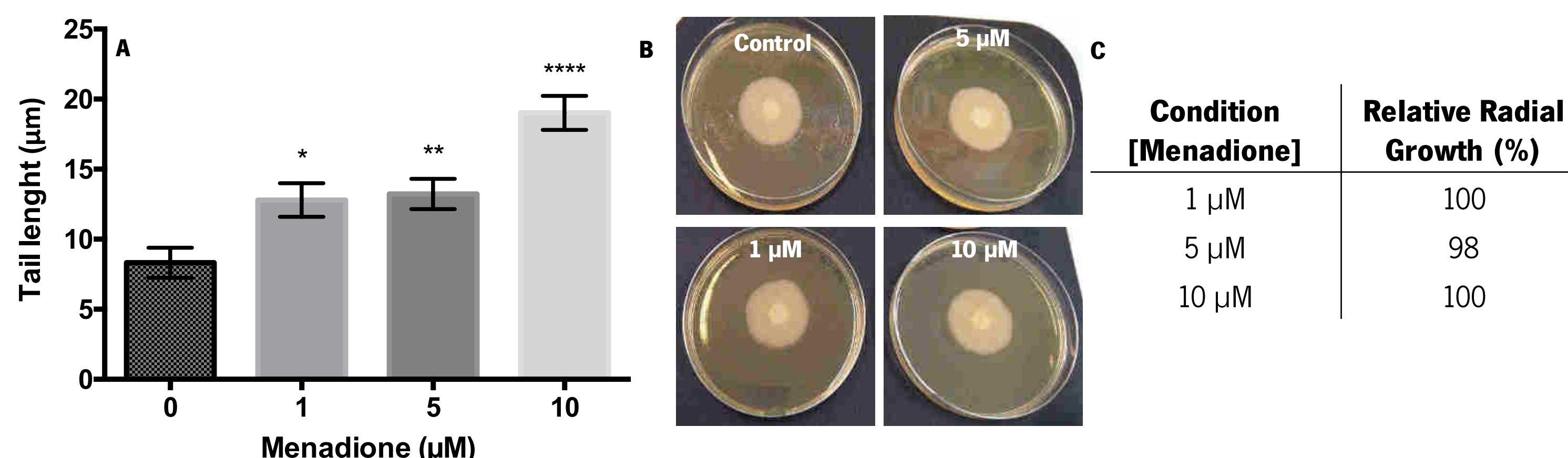


**Figure 2.2** – DNA damage and radial growth induced by exposure to Menadione. **(A)** DNA damage is represented as mean ± SEM of the tail length. **(B)** Representative photos of mycelia after growth for 3 days exposed to different concentrations of Menadione. **(C)** Relative growth of cells exposed to Menadione in comparison to the control condition (0 µM Menadione).

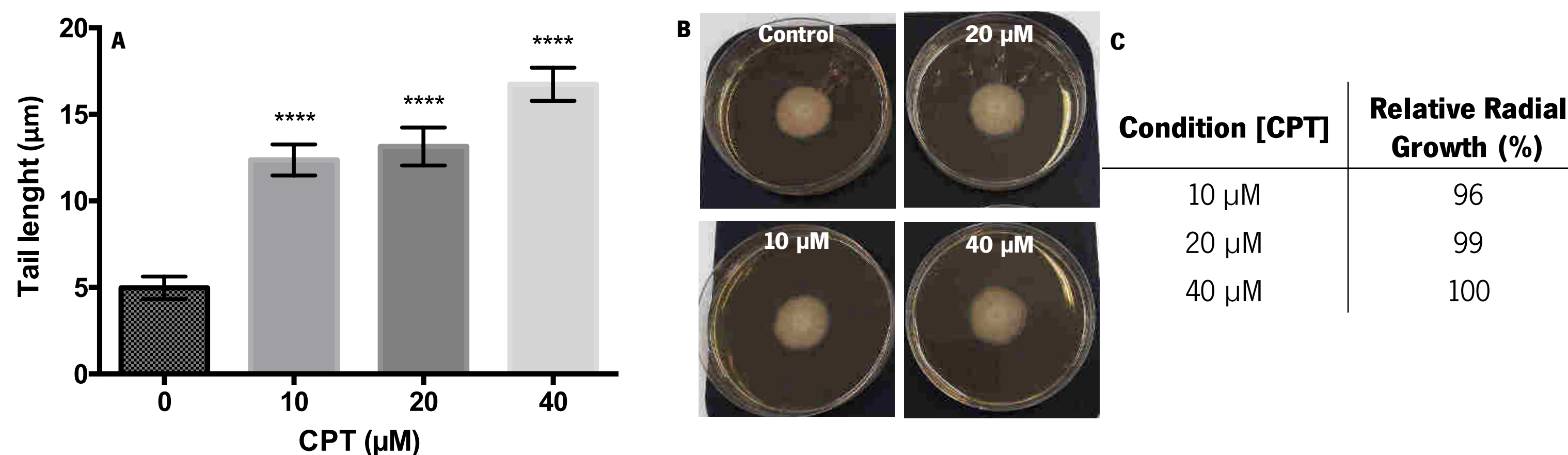

**Figure 2.3** – DNA damage and radial growth induced by exposure to CPT. **(A)** DNA damage is represented as mean ± SEM of the tail length. **(B)** Representative photos of mycelia after growth for 3 days exposed to different concentrations of CPT. **(C)** Relative growth of cells exposed to CPT in comparison to the control condition (0 µM CPT).

Statistic tests were determined by one-way ANOVA and Sidak's multiple comparisons test ( $p < 0.5$ ). (\*) represent statistically significant differences in comparison with the control condition (0). 3 biological replicates for each condition.

### 2. RIBOFLAVIN OVERPRODUCING STRAIN (A8) ONLY ACCUMULATES MORE DNA DAMAGE WHEN EXPOSED TO LIGHT DURING GROWTH

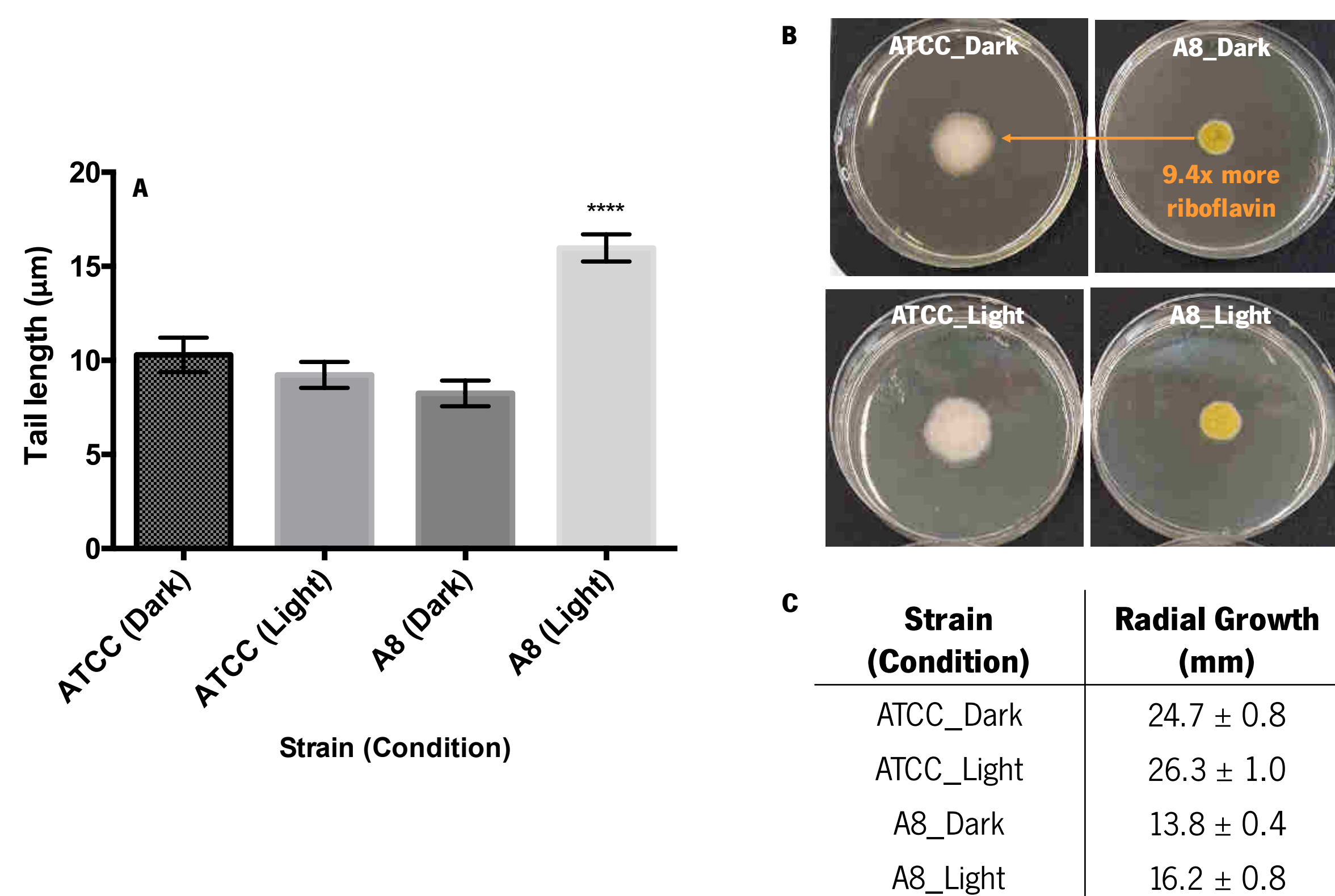

**Figure 3** – DNA damage and radial growth of different strains not exposed/exposed to light during the growth. **(A)** DNA damage is represented as mean ± SEM of the tail length. **(B)** Representative photos of mycelia after growth for 3 days in the different conditions tested. **(C)** Radial growth of mycelia in the different conditions tested (data are represented as mean ± SD).

## CONCLUSIONS

- We efficiently **adapted** the **Yeast Comet Assay** for use in *A. gossypii* with some alterations to the original protocol.
- This adapted protocol is a suitable tool to **measure oxidative and non-oxidative damage** in the **DNA** of *A. gossypii*, demonstrating **high sensitivity** and **response to different stress agents at different concentrations**.
- Using this method we report that **a riboflavin overproducing *A. gossypii* wild strain (A8)** accumulates **more DNA damage** than a non-overproducing strain when it is **exposed to light**, suggesting a **photosensitizer** role of **riboflavin** in these conditions, as already observed in mammalian cells<sup>4,5</sup>.

## REFERENCES

- Collins, Mol Biotechnol 2004; 26, 249-61.
- Langie et al., Front Genet 2015; 6:266
- Azevedo et al., Yeast 2011; 28, 55-61.
- Besaratinia, PNAS 2007; 104, 5953-58
- Alam, Arch Biochem Biophys 2015; 584, 10-9

## ACKNOWLEDGEMENTS

The authors thank the financial support from FCT, Portugal: strategic funding of UIDB/04469/2013 unit and COMPETE 2020 (POCI-01-0145-FEDER-005684) (Post-Doc fellowship to TQ Aguiar), BioTecNorte operation (NORTE-01-0145-FEDER-000004) funded by FEDER under the scope of Norte2020 - Programa Operacional Regional do Norte, and PhD grant to R Silva (PD/BD/113812/2015) through the FCT PhD Programme in Applied and Environmental Microbiology.