

Evaluating the Role of Host AMPK in *Leishmania* Burden 2

Diana Moreira, Jérôme Estaquier, Anabela Cordeiro-da-Silva, 3
and Ricardo Silvestre 4

Abstract 5

The study of host AMP-activated protein kinase (AMPK) activation during *Leishmania* infection imposes 6
distinct types of techniques to measure protein expression and activation, as well as to quantify, at 7
transcription and translational levels, its downstream targets. The investigation of host AMPK protein 8
modulation during *Leishmania* infection should primarily be assessed during in vitro infections using as a 9
host murine bone marrow-derived macrophages (BMMos). The infection outcome is assessed measuring 10
the percentage of infected cells in the context of BMMos. To evaluate AMPK activity during infection, the 11
expression of AMPK-phosphorylated at Thr172 as well as the transcription and translational levels of its 12
downstream targets are evaluated by quantitative PCR and immunoblotting. The modulation of AMPK 13
activity in vivo is determined specifically in sorted splenic macrophages harboring *Leishmania* parasites 14
recovered from infected mice using fluorescent-labeled parasites in the infectious inoculum. The modula- 15
tion of AMPK activity was assessed by AMPK activators and inhibitors and also using AMPK, SIRT1, or 16
LKB1-KO mice models. The infection outcome in BMMos and in vivo was further determined using these 17
two different approaches. To finally understand the metabolic impact of AMPK during infection, in vitro 18
metabolic assays in infected BMMos were measured in the bioenergetic profile using an extracellular flux 19
analyzer. 20

Key words *Leishmania*, AMPK, Bioenergetic profile, Extracellular flux analyzer, AMPK activators and 21
inhibitors, SIRT1, Mitochondria, Cell metabolism, Macrophages 22

1 Introduction 23

Leishmania spp. is the causative agent of leishmaniasis, a neglected 24
tropical disease transmitted by the bite of an infected female sand fly 25
[1]. These parasites are mainly phagocytosed by macrophages being 26
able to subvert their intracellular signaling pathways and compete 27
for similar resources [2–4]. A key question in the context of host- 28
pathogen interactions is how pathogens survive in a hostile envi- 29
ronment and how they hijack host machinery for their own benefit. 30

To address AMPK signaling during *Leishmania* infection, bone 31
marrow precursors are differentiated in vitro with macrophage 32
colony-stimulating factor (M-CSF) and used as a target cell for 33

infection (Subheading 3.1). This differentiation procedure leads to a higher yield and reproducibility, being these bone marrow-derived macrophages (BMMos) commonly considered as a model for the role of resident macrophages [5]. The infection is performed at early (6 h post-infection) and at later time points of infection (18 h post-infection) to acquire a dynamic profile of the infection process. The level of infected cells was obtained through the infection of BMMos with CFSE-labeled parasites being the percentage of CFSE⁺ or eFluor670⁺ cells, obtained by FACS analysis, a direct measurement of parasites internalization [6, 7] (Subheading 3.1). A direct analysis of AMPK activation in a context of infection is addressed through immunoblotting of total and phosphorylated Thr172 both in infected BMMos (Subheading 3.1) and in sorted infected splenic macrophages (Subheading 3.2). The modulation of AMPK activity during infection can be obtained establishing a pharmacological approach where an AMPK activator (AICAR) and/or inhibitor (compound C) can be used isolated or in combination, as was described by us and by other authors in different contexts of infection [7–9]. The establishment on an in vivo infection (Subheading 3.2) using myeloid-restricted (Mac)-AMPK, SIRT1, or LKB1 KO mice (Subheading 3.3) is imperative to evaluate the impact of AMPK for the infection outcome in *Leishmania*-parasitized organs. SIRT1 has been investigated in different contexts as a potential modulator of AMPK activation. On one hand SIRT1 protein has been described as an upstream activator of AMPK through LKB1 deacetylation and on the other hand has been considered a downstream target of AMPK, becoming activated by the increase levels of NAD⁺ [10–13]. The parasite load in *Leishmania*-parasitized organs is determined ultimately by real-time quantitative PCR (qRT-PCR) (Subheading 3.2) [7, 14]. The metabolic impact of AMPK during *Leishmania* infection can be finally addressed using the extracellular flux analyzer in infected BMMos. Host bioenergetic profile at real time is traced at basal conditions and in response to distinct pharmacological agents, which allow the quantification of metabolic parameters as extracellular acidification rate (ECAR), oxygen consumption rate (OCR), spare respiratory capacity (SRC), and glycolytic capacity (Subheading 3.3). Overall, with these techniques we could trace the activation of host AMPK network during *Leishmania* infection and the impact on parasite survival.

2 Materials

2.1 Animals and Parasites

1. Myeloid cell-specific (Mac)-Sirt1 KO mice, Mac-AMPK α 1 KO, Mac-LKB1 KO on C57BL/6 genetic background, and the respective littermate lox control (Lysozyme M-Cre^{+/+}

	Sirt1 ^{flox/flox} mice. All animals used in experiments are aged from 6 to 12 weeks.	78 79
	2. <i>L. infantum</i> (MHOM/MA/67/ITMAP-263) promastigotes.	80 81
2.2 Culture Reagents	1. Macrophage medium: DMEM, 4.5 g/L glucose, 20 mM HEPES, 10% (v/v) heat-inactivated fetal bovine serum (FBS), 2 mM L-glutamine, 100 U/ml penicillin, 100 mg/ml streptomycin.	82 83 84 85
	2. <i>Leishmania</i> medium: RPMI 1640, 20 mM HEPES, 10% (v/v) heat-inactivated FBS, 2 mM L-glutamine, 100 U/ml penicillin, 100 mg/ml streptomycin.	86 87 88
	3. Macrophage growth factor: 100 µg/ml M-CSF.	89
	4. XF medium: unbuffered DMEM, 4.5 g/L glucose, 2% (v/v) FBS, 2 mM L-glutamine, 100 U/ml penicillin, 100 mg/ml streptomycin.	90 91 92
	5. Carboxyfluorescein succinimidyl ester (CFSE, FITC) and eFluor670 (APC): 5 mM stock solution.	93 94
	6. AICAR (5-Aminoimidazole-4-carboxamide ribonucleotide).	95
	7. Compound C (also known as Dorsomorphin).	96
	8. SRT1720.	97 98
2.3 FACS Staining	1. DMEM-EDTA: DMEM, 50 mM EDTA.	99
	2. Phosphate buffer saline (PBS): 145 mM NaCl, 2.7 mM KCl, 1.5 mM KH ₂ PO ₄ , 8 mM Na ₂ HPO ₄ · 2H ₂ O, pH 7.	100 101
	3. PBS-FBS: PBS, 2% (v/v) FBS.	102
	4. Probes: 14.6 mM 2-[N-(7-nitrobenz-2-oxa-1,3-diazol-4-yl) amino]-2-deoxy-D-glucose (2-NBDG) and 1 mg/ml 7-Aminoactinomycin D (7-AAD).	103 104 105 106
2.4 Antibodies	1. FACS: anti-mouse monoclonal antibodies—anti-CD11b-PE (M1/70), anti-Ly6C-PerCP/Cy5.5 (HK1.4), anti-Ly6G-Pacific Blue (1A8), and anti-F4/80-PerCP-Cy5.5 (BM8).	107 108 109
	2. Cell magnetic separation: CD3e and CD19 microbeads.	110
	3. PBS-EDTA buffer: PBS, 0.5% (w/v) BSA (bovine serum albumin), 2 mM EDTA, pH 7.2.	111 112
	4. Immunoblot: total AMPKα (23A3), AMPKα phosphorylated at Thr172, total SIRT1 (H-300), total PGC1β (E-9), β-actin (C4), total PGC1α (4C1.3), and horseradish peroxidase-coupled secondary reagents.	113 114 115 116 117
2.5 Quantitative PCR Analysis	1. TRIzol reagent.	118
	2. Chloroform.	119
	3. Isopropanol.	120

t.1 **Table 1**
t.2 **List of primers used for real-time quantitative PCR (qRT-PCR)**

t.2	Genes	Forward sequence	Reverse sequence	Accession number
t.3	<i>Ppargc1a</i>	AGCCGTGACCACTGACAACGAG	GCTGCATGGTTCTGAGTGCTAAG	NM_008904
t.4	<i>Slc2a4</i>	ACATACCTGACAGGGCAAGG	CGCCCTTAGTTGGTCAGAAG	NM_009204
t.5	<i>RPS29</i>	CACCCAGCAGACAGACAAACTG	GCACTCATCGTAGCGTTCCA	NM_009093

4. Ethanol 75% (c/v) in H₂O.121
5. H₂O RNase free.122
6. RNeasy micro kit.123
7. QIAmp DNA micro kit.124
8. cDNA synthesis Kit.125
9. SYBR Green Supermix.126
10. Primers: R221 and R332 primers parasite load quantification127
and primers for real-time quantitative PCR listed in Table 1.128
11. DNazol.129

2.6 Immunoblot

1. Lysis buffer: 50 mM Tris-HCL, pH 7.4, 1% (v/v) Triton131
X-100, 150 mM NaCl, 10% (v/v) glycerol, 50 mM NaF,132
5 mM sodium pyrophosphate, 1 mM Na₃VO₄, 25 mM133
sodium-β-glycerophosphate, 1 mM DTT, 0.5 mM PMSF, pro-134
tease and phosphatase inhibitor cocktails.135
2. Dc protein assay with reagent S (Bio-Rad).136
3. SDS polyacrylamide gel 10% resolving gel buffer (4.5 mL):137
1.9 mL H₂O, 1.7 mL 30% acrylamide, 1.3 mL Tris-HCL pH138
8.8, 50 μL of 10% SDS, 50 μL of 10% APS, 2 μL of TEMED.139
Stacking gel (3 mL): 2.1 mL of H₂O, 500 μL of 30% acrylam-140
ide, 380 μL of 0.5 M Tris-HCL pH 6.8, 30 μL of 10% SDS,141
30 μL of 10% APS, 3 μL of TEMED.142
4. Transfer pack with nitrocellulose membrane.143
5. Western blot transfer system.144
6. Ponceau S.145
7. Block/diluent solutions: 10% (w/v) BSA, 0.05% (v/v)146
Tween 20.147
8. SuperSignal West Pico or West Dura chemiluminescent148
substrate.149
9. Striping solution I: 0.2 M glycine, 0.5 M NaCl, pH 2.8.150
10. Striping solution II: 0.5 M Acetic acid, 0.5 M NaCl, pH 2.5.151
- 152

2.7 Equipment

1. Gamma irradiator.	153
2. Flow cytometer.	154
3. Western blot transfer system.	155
4. Microplate reader.	156
5. Western blot detection system.	157
6. Protein expression quantification software.	158
7. Spectrophotometer.	159
8. Real-time PCR thermal cycler.	160
9. Magnetic cell separator system.	161
10. Extracellular flux analyzer.	162
	163

3 Methods

164

3.1 Infection of Bone Marrow-Derived Macrophages with *L. infantum*

3.1.1 Isolation and Culture of Mouse Bone Marrow-Derived Macrophages

1. Anesthetize mice through isoflurane inhalation and euthanize by cervical dislocation.	165
2. Isolate the femurs and tibias of each mouse with a sterile scalp and scissor from the hind legs (<i>see Note 1</i>).	166
3. Holding the femurs and the tibias with the help of a scissor, cut off with a scalpel each tip of the bones. Recover bone marrow precursors by flushing the bone marrow with 3–5 mL of ice-cold complete macrophage medium (per femur or tibia) with the help of a syringe and a 26 G needle.	167
4. Centrifuge bone marrow cells at $300 \times g$ for 10 min at 4 °C and resuspend in complete macrophage medium.	169
5. Plate the bone marrow precursor cells without counting in a 75 cm ² culture flask in a volume of 15 mL of complete macrophage medium. Incubate at 37 °C with 5% CO ₂ for 4–6 h.	170
6. Discard the adherent cells (differentiated macrophages from the stroma), and recover the bone marrow precursors present in the supernatant. Count with trypan blue and seed the bone marrow precursors in macrophage medium with 20 ng/ml of M-CSF at a proportion of 1×10^5 cells in 200 µL of medium in 96-well plates, 2×10^5 cells in 400 µL of medium in 24-well plates, and 1×10^6 cells in 2 mL of medium in 6-well plates.	171
7. Renew M-CSF growth factor at day 4 of culture in each well by adding 20 ng/ml of M-CSF. Macrophages acquired a definitive differentiation status at day 7 of culture, being defined as bone marrow-derived macrophages (BMMos) (<i>see Note 2</i>).	172
1. Start a 5 mL culture of <i>L. infantum</i> promastigotes at a concentration of 1×10^6 promastigotes/ml of <i>Leishmania</i> medium (<i>see Note 3</i>).	173
	174
	175
	176
	177
	178
	179
	180
	181
	182
	183
	184
	185
	186
	187
	188
	189
	190
	191
	192
	193

2. At day 5 of growth, wash *L. infantum* promastigotes twice with 5 mL of PBS and resuspend the parasites in 1 mL of PBS. Count the parasites by diluting 1/10 in 2% (v/v) glutaraldehyde that fix the parasites. Dilute the *Leishmania* culture to 6×10^7 promastigotes in 1 mL of PBS. 194
195
196
197
198
3. Label with 5 μ M CFSE for 10 min or 1 μ M eFluor670 for 5 min at 37 °C followed by 5 min incubation at 4 °C to stop the reaction (*see Note 4*). 199
200
201
4. Wash the parasites twice with 5 mL of PBS, and resuspend in 1 mL of *Leishmania* medium before proceeding to infections. 202
203
5. Count the parasites by diluting 1/10 in 2% (v/v) glutaraldehyde and infect 7-day differentiated BMMos with labeled and unlabeled *L. infantum* promastigotes at a 1:10 ratio. 204
205
206
6. Irradiate 6×10^7 promastigotes suspended in 1 mL of *Leishmania* culture medium at 3000 Gy. Perform control experiments with those irradiated parasites using the previous co-culture ratio. 207
208
209
210
7. After 4 h of incubation, remove the medium of each well. Wash the cells, at least twice, with macrophage culture medium pipetting up and down several times. 211
212
213
8. For each round of washing, observe under the microscope if the parasites are still present in the supernatant. Repeat the washing procedure until the complete removal of non-internalized parasites. 214
215
216
217
9. Detach BMMos after 6 and 18 h post-infection. Remove the macrophage medium and add an equivalent volume of DMEM-EDTA solution. Wait 5 min at room temperature, and recover the BMMos by pipetting up and down several times. 218
219
220
221
222
10. Centrifuge at $300 \times g$ for 10 min at room temperature. 223
11. Resuspend BMMos in 200 μ L of PBS-2%FBS solution, and incubate for 15 min at room temperature with 7-AAD at 1 μ g/ml. 224
225
226
12. Transfer the cells to the cytometer tubes. Acquire the samples in a flow cytometer. 227
228
13. In the cytometer adopt the following gate strategy: exclude the death cells that are stained positively to 7-AAD. In the viable population (7-AAD⁻), gate the cells that stained positively to eFluor670 (eFluor670⁺) or CFSE (CFSE⁺). 229
230
231
232
14. Obtain the percentage of infected cells by the % of CFSE⁺ or eFluor670⁺ cells. Determine the cell viability by the percentage of negative 7-AAD stained cells (*see Note 5*). 233
234
235
236

3.2 Infection of Mice with *L. infantum* and Sorting of Infected Splenic Macrophages

3.2.1 In Vivo Infection

3.2.2 Parasite Load Quantification

3.2.3 Sorting of Infected Splenic Macrophages

1. Infect mice intraperitoneally, using a 26 G needle, with 1×10^8 CFSE-labeled *L. infantum* promastigotes resuspended in 200 μ L of sterile PBS.
1. At 10 days post-infection, anesthetize mice through isoflurane inhalation and euthanize by cervical dislocation. Remove and weigh spleen and liver. Transfer both organs to a 70 μ m mesh cell strainer, within a small petri dish, and use the syringe plunger to process the organs to generate a single cell suspension. Centrifuge at $300 \times g$ for 10 min at 4 $^{\circ}$ C, and resuspend in 5 mL of complete macrophage medium.
2. Extract DNA from 10 mg of spleen and liver (single cell suspensions) or 3×10^6 bone marrow cells by DNazol, according to the manufacturer instructions. Dissolve DNA in 100 μ L of nuclease-free water. Quantify the total DNA in a NanoDrop spectrophotometer and prepare a twofold serial concentrations dilution adjusted for each tissue.
3. Quantify *Leishmania infantum* DNA by qPCR using 1000 nM of R223 and 500 nM of R333 primers for the small subunit rRNA (SSUrRNA) [15]. As a template use 400 ng of total DNA in 20 μ L of reaction with SYBR Green Supermix, according to the manufacturer's instructions.
4. Perform a touchdown qPCR in Bio-Rad My Cycler iQ5, with a final annealing temperature of 65 $^{\circ}$ C [16]. The denaturation temperature is at 94 $^{\circ}$ C (5 s) and synthesis at 72 $^{\circ}$ C (10 s) with 30 cycles (see Note 6).
5. Extrapolate CTs from a standard curve constructed previously with a serial dilutions of *L. infantum* DNA (strain MHOM/MA/67/ITMAP-263) diluted in host DNA (from spleen of naïve mice). Calculate then *Leishmania* content expressed by parasites/ μ g of DNA (see Note 7).
1. Euthanize naïve and CFSE-*L. infantum* promastigote infected mice at 18 or 48 h post-infection. Collect the spleens into a falcon with 5 mL of macrophage culture medium.
2. Process the organs as **step 2** in Subheading 3.2.1. Determine the cell number.
3. Wash the cells twice with 5 mL PBS. Pipette off supernatant completely and resuspend cell pellet in 80 μ L of PBS-EDTA buffer per 10^7 cells.
4. Deplete splenic T and B lymphocytes using CD3 ϵ and CD19 microbeads coupled with depletion columns using a magnetic cell isolation separator.

5. Label the remaining cell suspension with 2 $\mu\text{g}/\text{ml}$ anti-CD11b-PE, anti-Ly6C-PerCP/Cy5.5 and anti-Ly6G-Pacific Blue diluted in 25 μL of PBS-FBS solution for 30 min at 4 $^{\circ}\text{C}$ in the dark. 280-283
6. Wash the cells twice with 200 μL of PBS-FBS solution, and resuspend the cell pellet in 200 μL of PBS-FBS solution. Transfer the cell suspension to the cytometer tubes. 284-286
7. Sort the cells according to the surface expression of CD11b⁺Ly6C^{int/high}Ly6G^{low} and CFSE expression gated on infected (CFSE⁺CD11b⁺Ly6C^{int/high}Ly6G^{low}) or bystander (CFSE⁻CD11b⁺Ly6C^{int/high}Ly6G^{low}) splenic macrophages. As a control, sort CD11b⁺Ly6C^{int/high}Ly6G^{low} cells from the spleen of non-infected mice. The purity of the separation should be higher than 90%. 287-293
8. Count the cells obtained after the sorting assay using trypan blue to exclude dead cells. 294-295
9. Sorted splenic macrophages are rinse with PBS by centrifugation at $300 \times g$ for 10 min at 4 $^{\circ}\text{C}$. 296-297
10. Resuspend 1×10^5 sorted macrophages in 350 μL of RLT buffer with 7 μL of 2-mercaptoethanol, and isolate RNA according to RNeasy micro kit manufacturer's instructions. 298-300
11. Analyze the transcription levels of Ppargc1a (PGC-1 α) and Slc2a4 (GLUT4) by qPCR using RPS29 as housekeeping gene. 301-303
12. Resuspend 1×10^6 sorted macrophages in ice-cold lysis buffer for 30 min at 4 $^{\circ}\text{C}$ with shaking. The supernatant is recovered after centrifugation at $17,000 \times g$ during 20 min at 4 $^{\circ}\text{C}$. 304-306
13. Analyze AMPK activation by immunoblot using the experimental procedure described in Chapter 27 (*see Note 8*). 307-309

3.3 Modulation of AMPK in a Context of Host-*L. infantum* interaction

1. Infect KO mouse models (Mac-Sirt1 KO mice, Mac-AMPK α 1 KO, Mac-LKB1) and the respective littermate lox controls, as described in Subheading 3.2. 310-312
2. Evaluate the parasite load by qPCR as described in Subheading 3.2. 313-314

3.3.1 Modulation of AMPK Activity During In Vivo *L. infantum* Infection

3.3.2 Modulation of AMPK Activity in In Vitro *L. infantum*-Infected BMMos

1. Treat BMMos, from WT, Mac-SIRT1 KO, Mac-AMPK KO, and Mac-LKB1 KO mice, previously infected with CFSE or eFluor670-*L. infantum* promastigotes, at 6 h post-infection with 440 μM AICAR, 440 μM AICAR, and 5 μM compound C, 1 μM SRT1720 or left untreated (*see Note 9*). 316-320

3.3.3 In Vitro Metabolic Assays of Infected Bone Marrow Macrophages

2. After 24 h of infection, detach BMMOs by using the DMEM-EDTA solution and centrifuge at $300 \times g$ for 10 min at room temperature.
3. Resuspend cells in 200 μ L of PBS-2%FBS solution and incubate for 15 min at room temperature with 7-AAD at 1 μ g/ml.
4. Evaluate the infection rate as described in the Subheading 3.1.
1. Collect bone marrow precursors as described in the Subheading 3.1. Perform the 7-day differentiation process in 75 cm^2 culture flasks.
2. After 7 days in culture, scrap the cells and seed BMMo at 2×10^5 cells/well in 400 μ L of complete macrophage medium in XF-24 cell culture plates. Let the cells to adhere during an overnight period (*see* Note 10).
3. This procedure warrants a homogeneous distribution of the cells under the XF-24 plate surface, decreasing the variability among wells.
4. After an overnight period, infect the cells with irradiated or not *L. infantum* promastigotes at a 1:10 ratio, during 6 and 18 h.
5. One hour before the defined times of infection, wash the cells with pre-warmed XF medium. In the final wash, add 200 μ L of XF medium to each well, and incubate for an hour at 37 $^\circ\text{C}$ without CO_2 .
6. Determine the real-time measurement of bioenergetic profile (eight wells per condition), oxygen consumption rate (OCR), and extracellular acidification rate (ECAR), under basal conditions and in response to oligomycin (1 μM), fluoro-carbonyl cyanide phenylhydrazone (FCCP—1 μM), rotenone (1 μM), and antimycin A (1 μM), using an extracellular flux analyzer, at 6 and 18 h post-infection (Fig. 1).

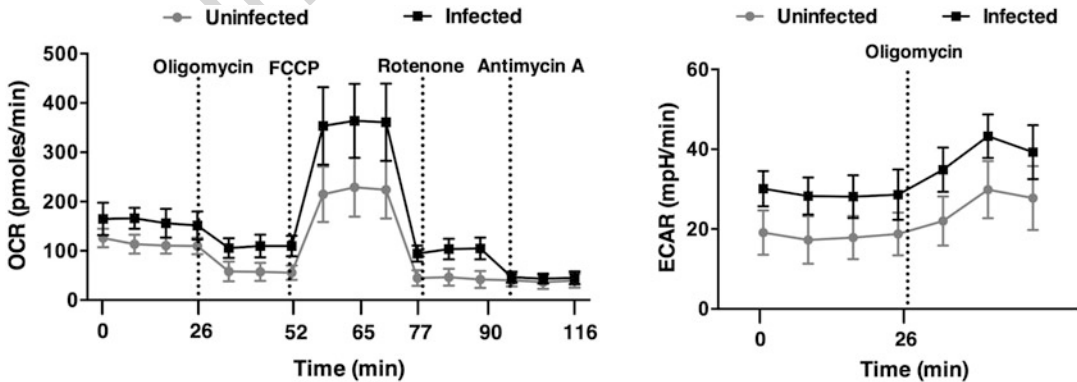


Fig. 1 A representative image of the real-time measurement of oxygen consumption rate (OCR) and the extracellular acidification rate (ECAR) in uninfected and infected BMMo. The different bioenergetic profiles are obtained under basal conditions and in response to oligomycin (1 μM), fluoro-carbonyl cyanide phenylhydrazone (FCCP—1 μM), rotenone (1 μM), and antimycin A (1 μM), using an XF-24 Extracellular Flux Analyzer

7. Obtain the non-mitochondrial respiration by subtracting the rotenone/antimycin A values. Calculate the spare respiratory capacity (SRC) by subtracting FCCP from basal OCR values, and define the glycolytic capacity as the variation between oligomycin and basal ECAR values. Normalize the values for all the conditions in relation to the protein content.

4 Notes

1. To maintain aseptic conditions and to easily remove the surrounding tissue, collect the femurs and tibias into alcohol solution for few seconds transferring afterward to DMEM medium.
2. Determine the purity of the macrophage culture by quantifying by flow cytometry the percentage of the cells expressing simultaneously the two surface markers CD11b and F4/80, which should be superior to 90%.
3. Maintain a cloned line of virulent *L. infantum* (MHOM/MA/67/ITMAP-263) by weekly subpassages at 26 °C in *Leishmania* medium. Use *L. infantum* promastigotes under four to ten passages in all the experiments.
4. eFluor 670 has a peak emission at 670 nm being detected with a 660/20 band-pass filter (equivalent to APC, Alexa Fluor™ 647, or eFluor™ 660) while CFSE peak emission of 521 nm being detected with a 530/30 nm band-pass emission filter (equivalent to GFP or FITC). Given that both dyes display similar results, the choice of labeling dye relates to the availability of the flow cytometry detectors and possibility to multi-color flow cytometry.
5. Giemsa staining is an alternative procedure to determine the percentage of infected cells. Upon **step 8**, Subheading **3.1.2**, remove the macrophage medium and wash BMMos with 1 mL of PBS. Add 1 mL of 1% (v/v) of paraformaldehyde for 20 min. Rinse twice with 1 mL of PBS and add 1 mL of Giemsa stain previously diluted 1:20 in deionized water. Incubate for 30 min and rinse with deionized water. Air dry and count the infected cells and the number of intracellular parasites in a microscope at a 400× magnification. Cell nucleus and the parasites acquire a purple coloration. The volumes given consider an infection protocol for 2×10^5 BMMos in 400 µL of medium in 24-well plates. Up or downscale accordingly to the chosen infection protocol.
6. Whenever the qPCR gave a positive (with the expected melting curve) but unquantifiable value or a doubtful specific product (aberrant melting curve), perform a nested PCR [17] that has a

- higher sensitivity (0.01 parasites) than the qPCR (0.6 parasites) to confirm the positivity of the quantitative result. For the first amplification reaction, use 300 nM of R221 and R332 primers [17]. For the second reaction, use 10 μ L of the first PCR product diluted 1:40, which will serve as a template with the same R223 and R333 primers (300 nM and 150 nM, respectively) used for the qPCR.
7. As alternatives, the detection and quantification of *L. infantum* kinetoplastid DNA can be performed by TaqMan-based qPCR assay [18] or the parasite burden determined by limiting dilution assay [19]. For the former, reaction mixtures are composed of ABI TaqMan PCR 2 \times (Applied Biosystems), 375 nM of direct primer (CTTTTCTGGTCCTCCGGGTAGG), 375 nM of reverse primer (CCACCCGGCCCTATTTTACACCAA), 250 nM of hydrolysis probe (5'FAM-TTTTCGCAGAACGCCCCTACCCGC-3'TAMRA), and 100 ng of sample DNA. Thermocycling settings consist of one hold of 10 min at 95 $^{\circ}$ C followed by a two-step temperature (95 $^{\circ}$ C for 15 s and 60 $^{\circ}$ C for 60 s) over 40 cycles. A standard curve is established corresponding to a range of 50,000–0.01 parasites. Sample normalization is performed by quantifying a host gene (murine albumin), in 10 μ L parallel reactions consisting of SYBR Green ROX Mix 2 \times , 100 nM of forward primer (CCATTGGTGAGACCAGAGGT), 100 nM of reverse primer (GAGGCAGGCAGCTTTATCAG), 100 ng of DNA, and the same thermal profile used for parasite quantification. A calibration curve ranging from 10,000 to 0.1 cells is established and parasite load expressed as the number of parasites per million of host cells. For the parasite burden, remove the organs from the mice (spleen and liver), weight and homogenize in RPMI medium. After cell counting, perform a subsequent twofold dilutions, in quadruplicate, in 96-well plates, and then incubate at 26 $^{\circ}$ C for 15 days. Record the presence or absence of motile promastigotes in each well. Calculate the number of parasites per gram of organ (parasite burden) as follows: parasite burden = [(geometric mean of reciprocal titer from each quadruplicate cell culture/weight of homogenized organ) reciprocal fraction of the homogenized organ inoculated into the first well].
 8. AICAR is an AMPK activator, compound c is an AMPK inhibitor, and SRT1720 is an activator of SIRT1. BMMos from Mac-SIRT1 KO mice are treated with AICAR and AICAR + compound c, and BMMos from Mac-AMPK KO and Mac-LKB1 KO mice are treated with SRT1720, respectively.
 9. In the immunoblot assays, use as a readout the expression of total AMPK, AMPK-PThr172, and PGC-1 α .

10. The procedure warrants a homogeneous distribution of the cells under the XF-24 plate surface, decreasing the variability among wells.

Acknowledgments

This work was supported by the Northern Portugal Regional Operational Programme (NORTE 2020), under the Portugal 2020 Partnership Agreement, through the European Regional Development Fund (FEDER) (NORTE-01-0145-FEDER-000013) and the Fundação para a Ciência e Tecnologia (FCT) (contract IF/00021/2014 to R.S.), by FEDER funds through the Operational Competitiveness Programme (COMPETE), and by national funds through FCT (Fundação para a Ciência e a Tecnologia) under the project FCOMP-01-0124-FEDER-011054 (PTDC/SAU-FCF/100749/2008) and PTDC/BIA-MIC/118644/2010. The research leading to these results has also received funding from the European community's Seventh Framework Programme under grant agreement No.602773 to JE and ACS (Project KINDRED). DM was supported by SFRH/BD/91543/2012. JE thanks the Canada Research Chair program for their supports.

References

1. Murray HW, Berman JD, Davies CR, Saravia NG (2005) Advances in leishmaniasis. *Lancet* 366:1561–1577
2. Estaquier J, Vallette F, Vayssiere JL, Mignotte B (2012) The mitochondrial pathways of apoptosis. *Adv Exp Med Biol* 942:157–183
3. Rodrigues V, Cordeiro-da-Silva A, Laforge M, Ouaisi A, Silvestre R, Estaquier J (2012) Modulation of mammalian apoptotic pathways by intracellular protozoan parasites. *Cell Microbiol* 14:325–333
4. Rodrigues V, Cordeiro-da-Silva A, Laforge M, Ouaisi A, Akharid K, Silvestre R, Estaquier J (2014) Impairment of T cell function in parasitic infections. *PLoS Negl Trop Dis* 8:e2567
5. Manzanero S (2012) Generation of mouse bone marrow-derived macrophages. *Methods Mol Biol* 844:177–181
6. Resende M, Moreira D, Augusto J, Cunha J, Neves B, Cruz MT, Estaquier J, Cordeiro-da-Silva A, Silvestre R (2013) *Leishmania*-infected MHC class II high dendritic cells polarize CD4+ T cells toward a nonprotective T-bet+ IFN-gamma+ IL-10+ phenotype. *J Immunol* 191:262–273
7. Moreira D, Rodrigues V, Abengozar M, Rivas L, Rial E, Laforge M, Li X, Foretz M, Viollet B, Estaquier J, Cordeiro da Silva A, Silvestre R (2015) *Leishmania infantum* modulates host macrophage mitochondrial metabolism by hijacking the SIRT1-AMPK axis. *PLoS Pathog* 11:e1004684
8. Bulusu V, Thakur SS, Venkatachala R, Balaram H (2011) Mechanism of growth inhibition of intraerythrocytic stages of *Plasmodium falciparum* by 5-aminoimidazole-4-carboxamide ribonucleoside (AICAR). *Mol Biochem Parasitol* 177:1–11
9. Kondratowicz AS, Hunt CL, Davey RA, Cherry S, Maury WJ (2013) AMP-activated protein kinase is required for the macropinocytic internalization of ebolavirus. *J Virol* 87:746–755
10. Fulco M, Cen Y, Zhao P, Hoffman EP, McBurney MW, Sauve AA, Sartorelli V (2008) Glucose restriction inhibits skeletal myoblast differentiation by activating SIRT1 through AMPK-mediated regulation of Nampt. *Dev Cell* 14:661–673
11. Hou X, Xu S, Maitland-Toolan KA, Sato K, Jiang B, Ido Y, Lan F, Walsh K, Wierzbicki M,

- 515 Verbeuren TJ, Cohen RA, Zang M (2008) 541
516 SIRT1 regulates hepatocyte lipid metabolism 542
517 through activating AMP-activated protein 543
518 kinase. J Biol Chem 283:20015–20026 544
- 519 12. Lan F, Cacicedo JM, Ruderman N, Ido Y 545
520 (2008) SIRT1 modulation of the acetylation 546
521 status, cytosolic localization, and activity of 547
522 LKB1. Possible role in AMP-activated protein 548
523 kinase activation. J Biol Chem 549
524 283:27628–27635 550
- 525 13. Canto C, Gerhart-Hines Z, Feige JN, 551
526 Lagouge M, Noriega L, Milne JC, Elliott PJ, 552
527 Puigserver P, Auwerx J (2009) AMPK regu- 553
528 lates energy expenditure by modulating 554
529 NAD⁺ metabolism and SIRT1 activity. Nature 555
530 458:1056–1060 556
- 531 14. Cunha J, Carrillo E, Sanchez C, Cruz I, 557
532 Moreno J, Cordeiro-da-Silva A (2013) Char- 558
533 acterization of the biology and infectivity of 559
534 *Leishmania infantum* viscerotropic and dermo- 560
535 tropic strains isolated from HIV⁺ and HIV- 561
536 patients in the murine model of visceral leish- 562
537 maniasis. Parasit Vectors 6:122 563
- 538 15. van Eys GJ, Schoone GJ, Kroon NC, Ebeling 564
539 SB (1992) Sequence analysis of small subunit 565
540 ribosomal RNA genes and its use for detection 566
and identification of *Leishmania* parasites. Mol 541
Biochem Parasitol 51:133–142 542
16. Miro G, Oliva G, Cruz I, Canavate C, 543
Mortarino M, Vischer C, Bianciardi P (2009) 544
Multicentric, controlled clinical study to evalu- 545
ate effectiveness and safety of miltefosine and 546
allopurinol for canine leishmaniasis. Vet Der- 547
matol 20:397–404 548
17. Cruz I, Chicharro C, Nieto J, Bailo B, 549
Canavate C, Figueras MC, Alvar J (2006) 550
Comparison of new diagnostic tools for man- 551
agement of pediatric Mediterranean visceral 552
leishmaniasis. J Clin Microbiol 44:2343–2347 553
18. Rodrigues V, Laforge M, Campillo-Gimenez L, 554
Soundaramourty C, Correia-de-Oliveira A, 555
Dinis-Oliveira RJ, Ouaiissi A, Cordeiro-da- 556
Silva A, Silvestre R, Estaquier J (2014) Abort- 557
ive T follicular helper development is asso- 558
ciated with a defective humoral response in 559
Leishmania infantum-infected macaques. 560
PLoS Pathog 10:e1004096 561
19. Silvestre R, Cordeiro-Da-Silva A, Santarem N, 562
Vergnes B, Sereno D, Ouaiissi A (2007) SIR2- 563
deficient *Leishmania infantum* induces a defined 564
IFN-gamma/IL-10 pattern that correlates with 565
protection. J Immunol 179:3161–3170 566